


ragmask v2

ragmask.com


1. Print just this page.

You don't need the other pages, save your paper.

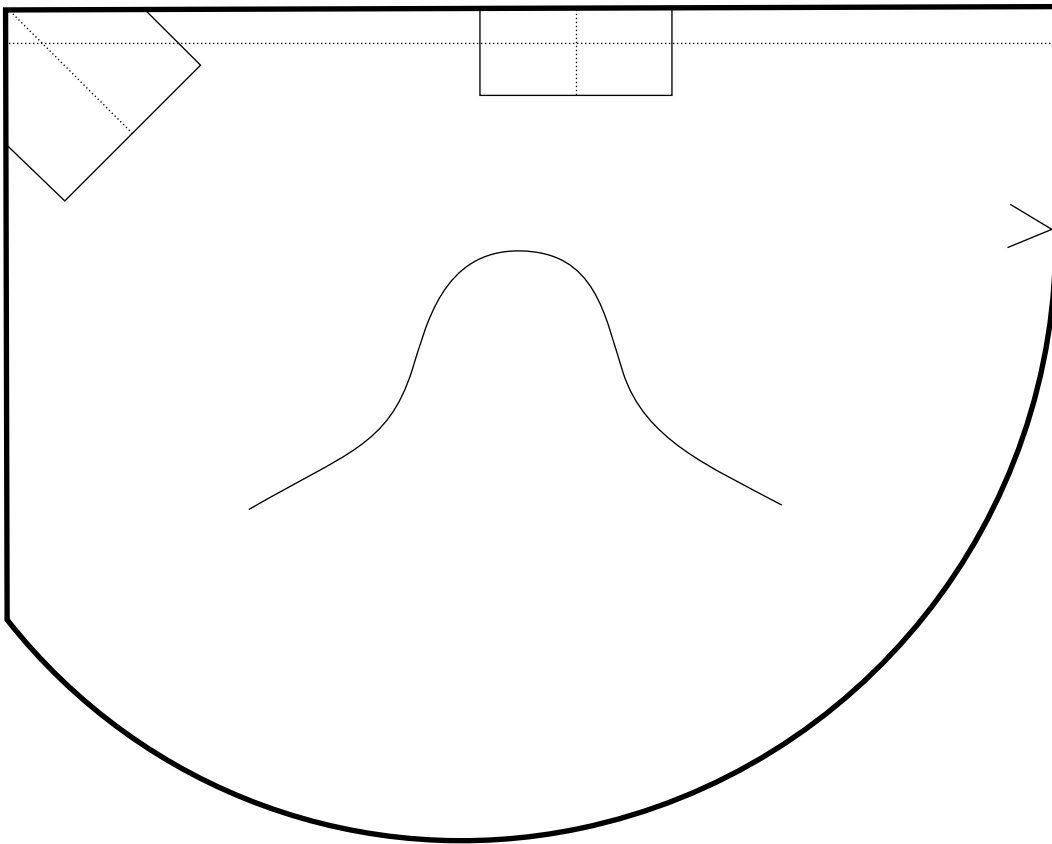

2. Check the scale.

Make sure the page printed at the right size, a real ruler should match the ruler below.

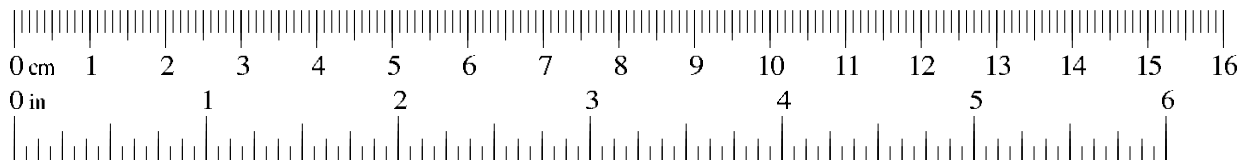

3. Cut materials.

- **11"x11"** (28x28cm) fabric, **1** piece per mask
- **1"x16"** (2.5x40cm) fabric for straps, **4** pieces per mask



- **10"x5"** (26x13cm) filter*, **1** piece per mask
- **4"** (10cm) metal wire, **1** piece per mask, 14ga copper wire works well

** effectiveness of various filter materials forthcoming...*


4. Fold fabric over filter.



5. Sew around 3 sides.
~3/16" (0.5cm) seam allowance, quick > pretty.




6. Fold.

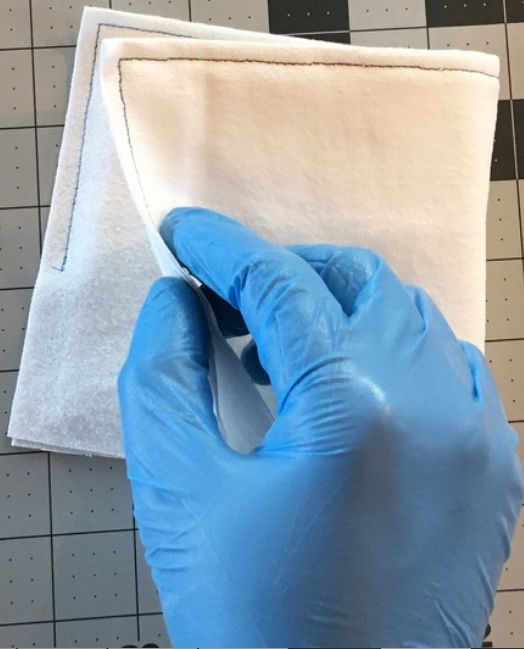

7. Mark.

Use template from page 1.

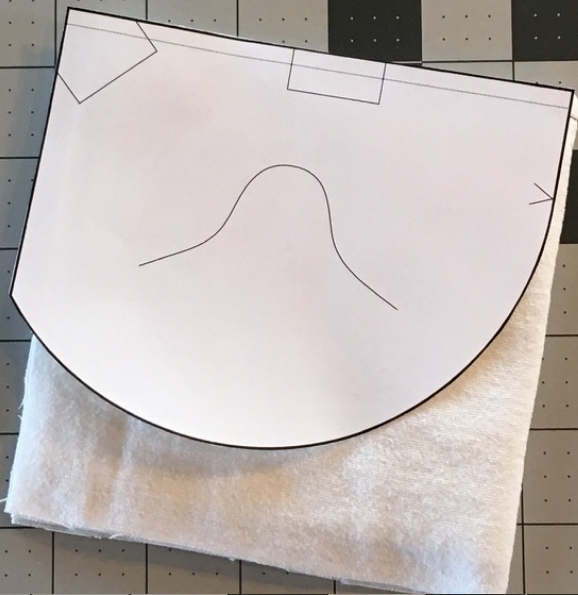

8. Sew.


9. Cut.



10. Flip inside out.




11. Mark strap locations.

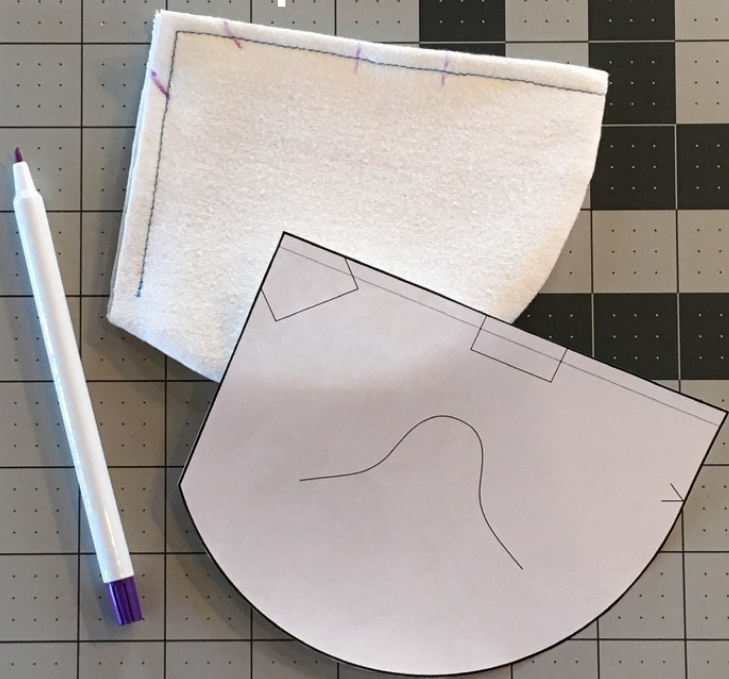


12. Sew straps on outside face of mask.



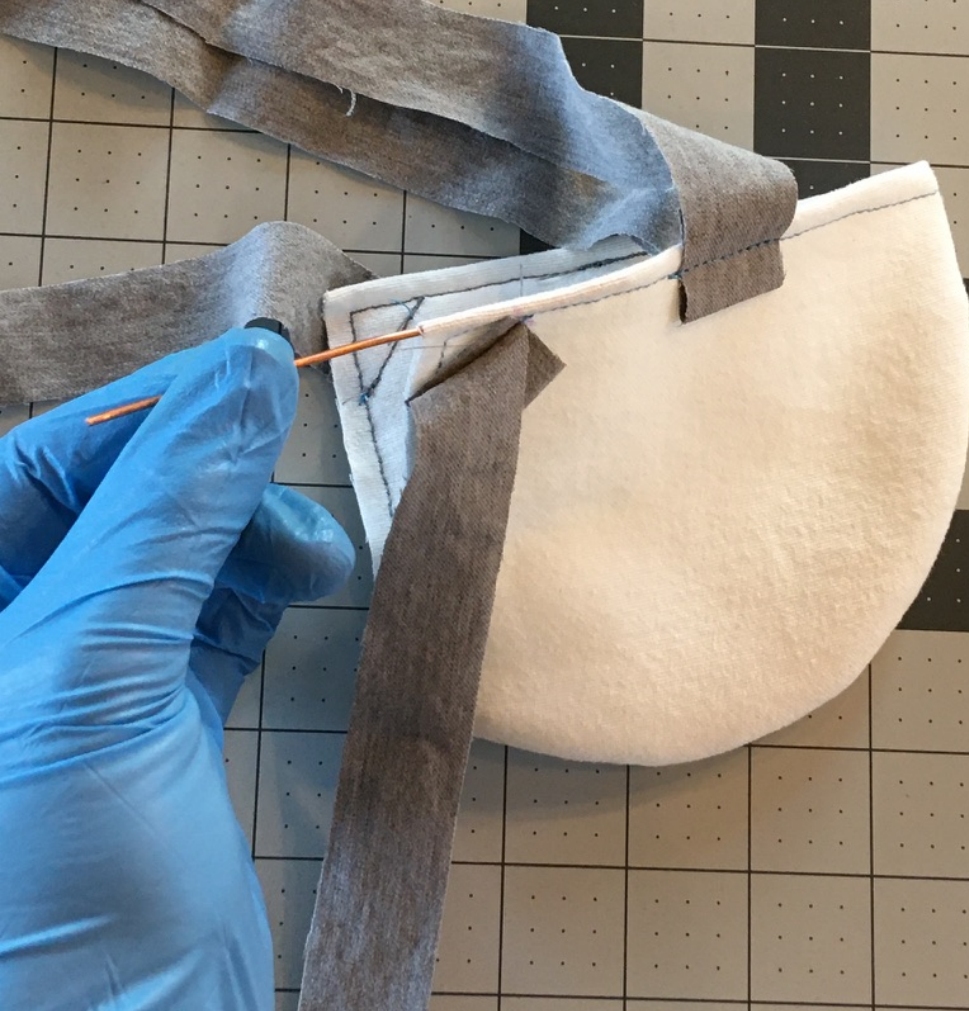

13. Insert wire and center over nose.

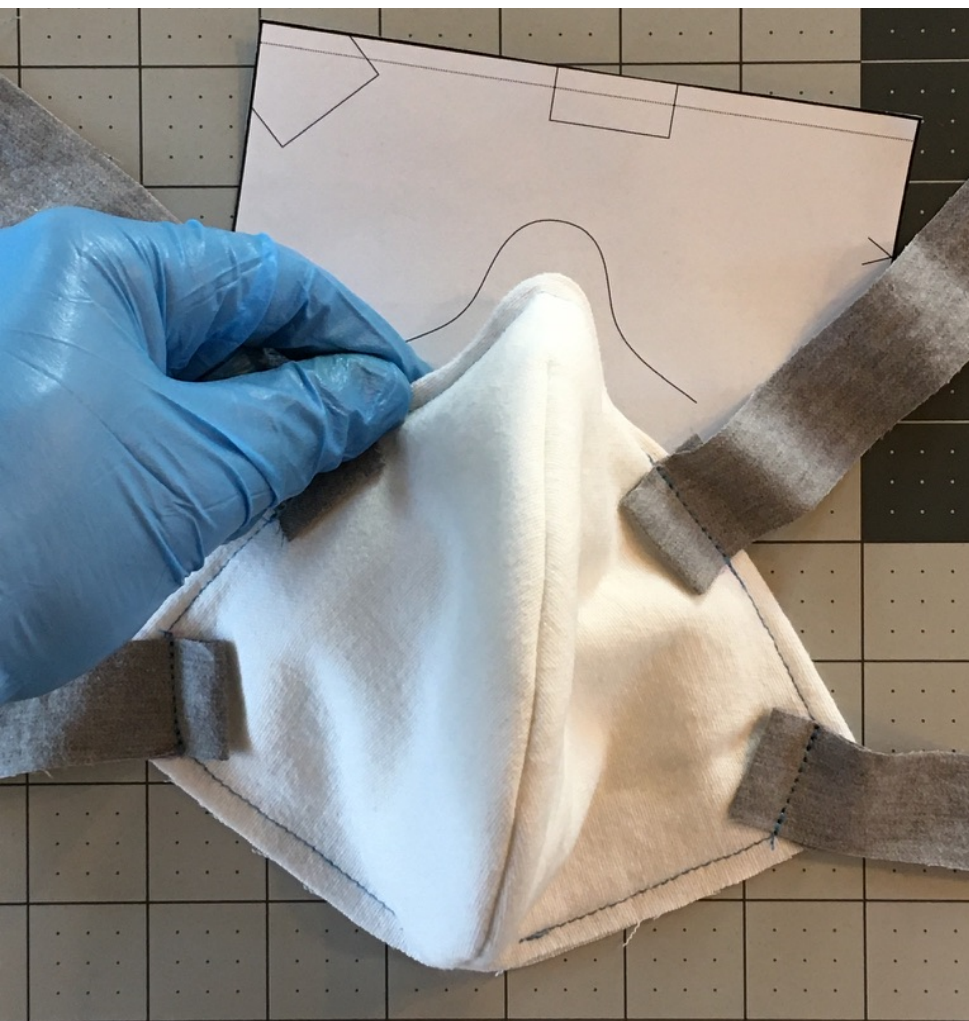

14. Form wire into approximate nose shape.

Wear

ragmask.com



Use and care details forthcoming...